©AMERIPHONE products are manufactured by Walker, a Division of Plantronics, Inc.
12082 Western Avenue, Garden Grove, CA 92841
(800) 874-3005 VOICE • (800) 772-2889 TTY/TDD • (714) 897-4703 FAX
email: ameriphonecs@plantronics.com website: www.ameriphone.com

The Extra Loud & Clear

XL-30 AMPLIFIED TELEPHONE

(Line Powered)

Operating Instructions

Important Safety Instructions

When using your telephone equipment, basic safety precautions should always be followed to reduce the risk of fire, electric shock and persons including the following:

1. Read and understand all instructions.
2. Follow all warnings and instructions marked on the telephone.
3. Do not use this telephone near a bathtub, wash basin, kitchen sink or laundry tub, in a wet basement, near a swimming pool or anywhere else there is water.
4. Avoid using a telephone (other than a cordless type) during a storm. There may be a remote risk of electrical shock from lightning.
5. Do not use the telephone to report a gas leak in the vicinity of the leak.
6. Unplug this telephone from the wall outlets before cleaning. Do not use liquid cleaners or aerosol cleaners on the telephone. Use a damp cloth for cleaning.
7. Place this telephone on a stable surface. Serious damage and/or injury may result if the telephone falls.
8. Do not cover the slots and openings on this telephone. This telephone should never be placed near or over a radiator or heat register. This telephone should not be placed in a built-in installation unless proper ventilation is provided.
9. Operate this telephone using the electrical voltage as stated on the base unit or the owner's manual. If you are not sure of the voltage in your home, consult your dealer or local power company.
10. Do not place anything on the power cord. Install the telephone where no one will step or trip on the cord.
11. Do not overload wall outlets or extension cords as this can increase the risk of fire or electrical shock.
12. Never push any objects through the slots in the telephone. They can touch dangerous voltage points or short out parts that could result in a risk of fire or electrical shock. Never spill liquid of any kind on the telephone.
13. To reduce the risk of electrical shock, do not take this phone

Important Safety Instructions

apart. Opening or removing covers may expose you to dangerous voltages or other risks. Incorrect reassembly can cause electric shock when the appliance is subsequently used.

14. Unplug this product from the wall outlets and refer servicing to the manufacturer under the following conditions:

- A. When the power supply cord or plug is frayed or damaged.
- B. If liquid has been spilled into the product.
- C. If the telephone has been exposed to rain or water.
- D. If the telephone does not operate normally by following the operating instructions. Adjust only those controls that are covered by the operating instructions. Improper adjustment may require extensive work by a qualified technician to restore the telephone to normal operation.
- E. If the telephone has been dropped or the case has been damaged.

- F. If the telephone exhibits a distinct change in

performance.

- 15. Never install telephone wiring during a lightning storm.
- 16. Never install telephone jacks in wet locations unless the jack is specifically designed for wet locations.
- 17. Never touch uninsulated telephone wires or terminals unless the telephone line has been disconnected at the network interface.
- 18. Use caution when installing or modifying telephone lines.
- 19. Use only the power cord and batteries indicated in this manual. Do not dispose of batteries in a fire. They may explode. Check with local codes for possible special disposal instructions.

ADDITIONAL SAFETY NOTES FOR CANADIAN USERS

The following items are included as part of the CS-03 Requirements. The standard connecting arrangement for the equipment is CA11A. This product meets the applicable Industry Canada technical specifications.

NOTICE: The Canadian

Important Safety Instructions

Department of Communications label identifies certified equipment. This certification means that the equipment meets certain telecommunications network protective operational and safety requirements. The Department does not guarantee that the equipment will operate to the user's satisfaction. Before installing this equipment, users should ensure that it is permissible to be connected to the facilities of the local telecommunications company. The equipment must also be installed using an acceptable method of connection. In some cases, the company's inside wiring associated with a single line individual service may be extended by means of a certified connector assembly (telephone extension cord). The customer should be aware that compliance with the above conditions may not prevent degradation of service in some situations. Repairs to certified equipment should be made by an authorized Canadian maintenance facility designated by the supplier. Any repairs or alterations made by the user to this equipment, or equipment malfunctions, may give

the telecommunications company cause to request the user disconnect the equipment.

Users should ensure for their own protection that the electrical ground connections of the power utility, telephone lines and internal metallic water pipe system, if present, are connected together. This precaution may be particularly important in rural areas.

CAUTION: Users should not attempt to make such connections themselves, but should contact the appropriate electric inspection authority, or electrician, as appropriate. The Ringer Equivalent Number is an indication of the maximum number of terminals allowed to be connected to a telephone interface. The termination on an interface may consist of any combination of devices subject only to the requirement that the sum of the Ringer Equivalent Number of all the devices not exceed five.

SAVE THESE INSTRUCTIONS

Introducing the XL-30

Warranty Service	1
Sales Receipt	1
Help from Ameriphone	1

Features	3
-----------------------	---

Setting up Your XL-30

Installing the Backup Batteries	4
Setting the Dialing Mode	4
Connecting for Desk or Wall Mount	4
Desktop Use	4
For Wall Mounting	5

Using Your XL-30

Incoming Voice Volume	6
Incoming Voice Tone	7
Ringer Volume and Pitch	8
Hold	8
Last Number Redial	8
Flash	8

Special Features

Programming the Memory Buttons	9
Direct Audio Output	10
Outgoing Voice Amplification	10
Cochlear Implant Adapter Cord	10

Troubleshooting Chart	13
------------------------------------	----

Regulatory Compliance	13
------------------------------------	----

Warranty	19
-----------------------	----

Specifications	21
-----------------------------	----

Index	22
--------------------	----

Introducing the XL-30

Thank you for selecting the XL-30 amplified telephone from Ameriphone. These Operating Instructions and the associated Quick Operating Guide provide you with the information you need to use your XL-30 effectively and safely. Read this manual thoroughly before using your telephone. Keep the manual near the telephone for easy reference.

Warranty Service

Your telephone is designed to provide years of quality service. But, should the phone malfunction and the Trouble-shooting Chart on page 11 not resolve the problem, follow the Warranty procedure on page 12.

Sales Receipt


Be sure to save your sales receipt as proof of purchase date should you need warranty service.

Help from Ameriphone


For help with using your XL-30, call our Customer Relations department at 800-874-3005.

Please make sure your XL-30 package includes the items shown in Figure 1.


Introducing the XL-30

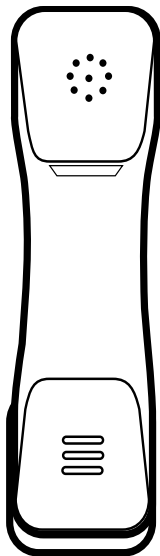Handset

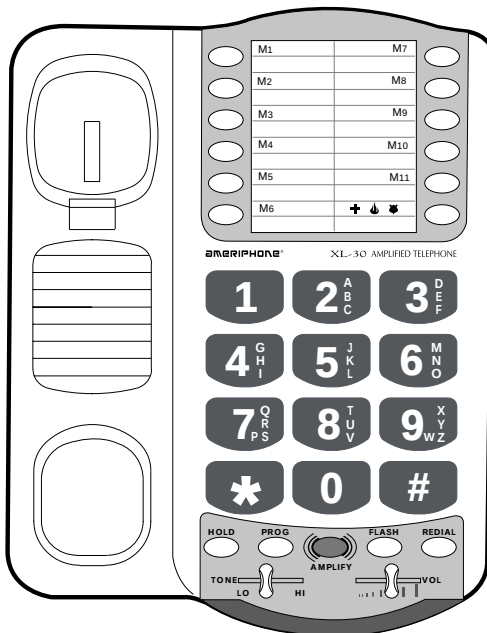
Base Unit

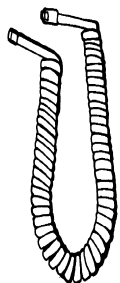
Handset
Coil Cord

Telephone Line Cords

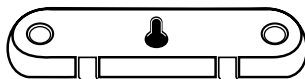
Mounting Bracket

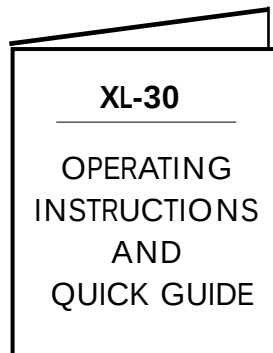Operating Guides

Figure 1 - Package Components

Features on the XL-30

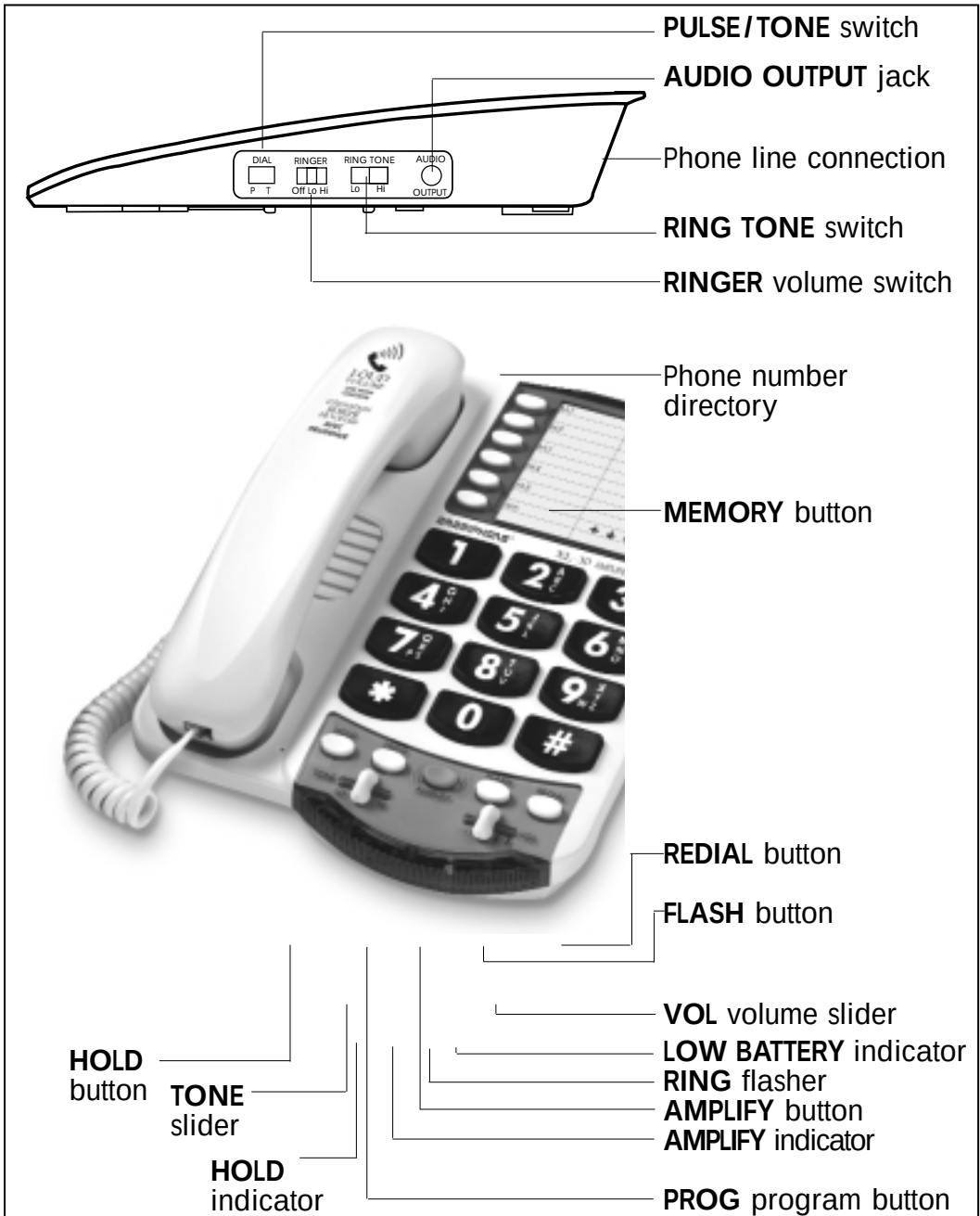3 Figure 2 - Base Unit Controls

Setting up Your XL-30

There are five steps involved in setting up your XL-30.

- Insert one 9 Volt alkaline battery for extra bright ring flasher (optional).
- Decide if you want the phone to sit on a desk or hang on the wall.
- Connect the telephone components.
- Program up to to 12 telephone numbers to call with the press of a memory button.
- Setup special features

Installing the Backup Batteries

To add extra brightness to the XL-30's ring flasher when you receive calls, install one 9 Volt alkaline battery (not included).

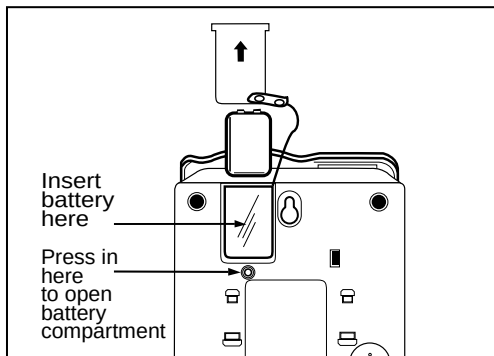
Figure 3 - Accessing the Batteries

To install new battery:

1. Slide open the battery compartment cover on the bottom of the telephone, (Figure 3).
2. Install one fresh battery, and reinstall the cover.

Setting the Dialing Mode

Set the dial mode switch to **T** (tone) if you have touch tone service. Set the switch to **P** (pulse) if you have rotary dialing.

Connecting for Desk or Wall Mount

Desktop Use

1. Connect the telephone line cord and the handset cord to the telephone as shown in Figure 4.
2. Lift the handset and listen for a dial tone.

The phone is ready to use.

Setting up Your XL-30 (continued)

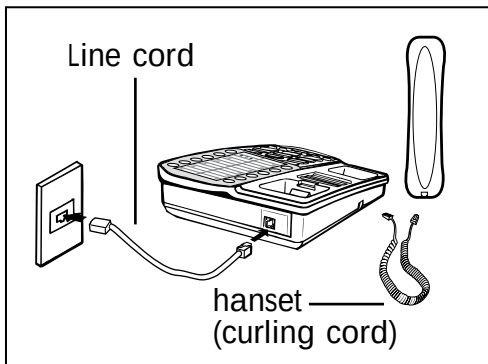
Figure 4 - Connecting the Components

For Wall Mounting

1. Pull the handset hook out as shown in Figure 5 and reverse its position so that the hook points up and will hold the handset when you hang up.

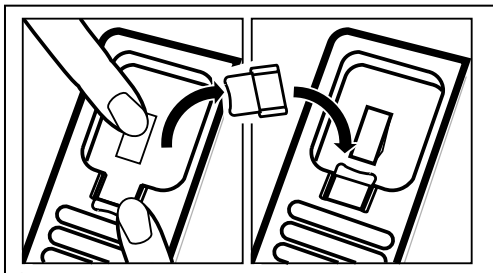
Figure 5 - Reversing the Hook

2. Insert the tabs of the phone base into the slots in the

lower position of the base as shown in Figure 6. Push the mounting bracket in and down firmly until it snaps into place.

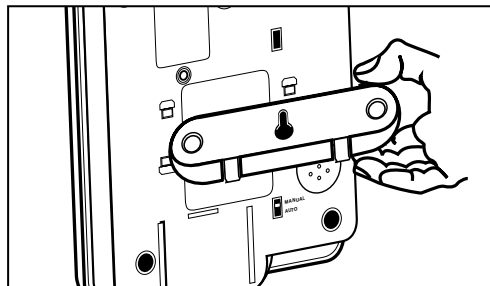
Figure 6 - Mounting Bracket for Wall Mounting

3. Locate a desired wall phone jack and hold phone upright.
4. Plug in the short phone line cord and mount the phone as shown in Figure 7.

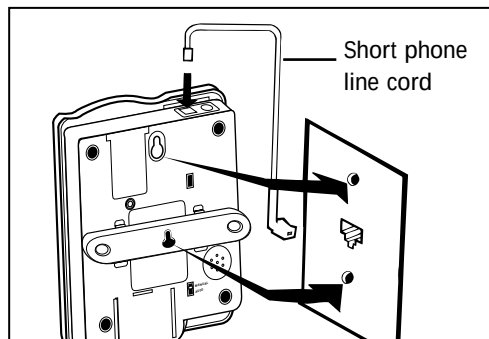
Figure 7 - Wall mounting installation

Setting up Your XL-30 (continued)

5. Angle the phone downward to feel for the LOWER protruding nail head. Insert the WALL JACK'S nail head into the LOWER part of the phone's mounting bracket.
6. Once the lower nail head has been inserted, insert the UPPER nail head into the phone's bracket and slide the phone down until it is firmly in place on the wall.
7. Once firmly in place, attach the short line cord to the phone and attach the phone's handset.
8. Lift the handset and listen for a dial tone. The phone is now ready to use.

Using Your XL-30

Incoming Voice Volume

You can adjust the volume of incoming calls by simply sliding the **VOL** control. This gives you up to 18 dB. more volume. If you want more amplification, follow the directions below.

1. Press the **AMPLIFY** button as shown in Figure 8 to turn the incoming voice amplifier ON or OFF. When **AMPLIFY** is on, the **AMPLIFY** indicator light comes on.

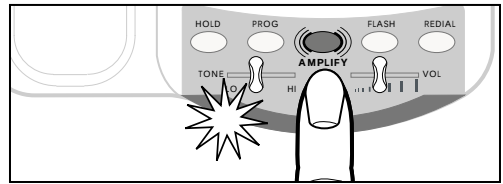
Figure 8 - AMPLIFY Button and Indicator

2. Adjust the volume in the ear piece by moving the **VOL** slide piece by moving the **VOL** slide control in Figure 9. With **AMPLIFY** ON incoming volume will be up to 30dB louder.

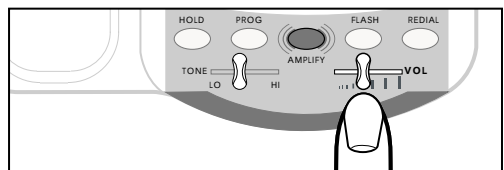
Figure 9 - Volume Slide Control

CAUTION: Repeated incremental exposure to amplification levels greater than 18 dB may be harmful to individuals without hearing disabilities. Therefore, **do not remove the warning label attached to the back of the handset.** If it is likely that a visually impaired person will use the telephone, securely attach the provided warning printed in Braille to the back of the handset (See, Diagram A).

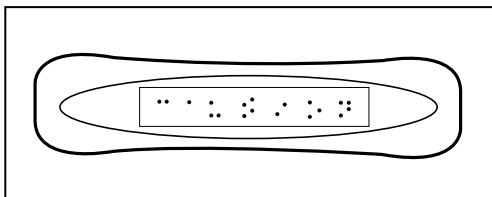
Diagram A.

NOTE: In standard use, the amplifier turns off whenever you hang up the telephone. This is a useful feature if many people use the phone. If you want the amplifier to always remain on, slide the Volume Reset Override switch on the back of the phone to **ON**, as shown in Figure 10.

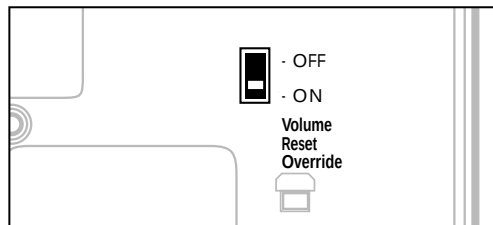
Figure 10 - Volume Reset Override On/Off

Incoming Voice Tone

The XL-30 provides extra amplification at the sound frequency you need to boost. To identify the sound frequency range that best suits your needs, follow these steps:

1. When you hear a voice on the line, press the **AMPLIFY** button shown in Figure 8.
2. Adjust the **TONE** slide control in Figure 11 to suit your hearing.

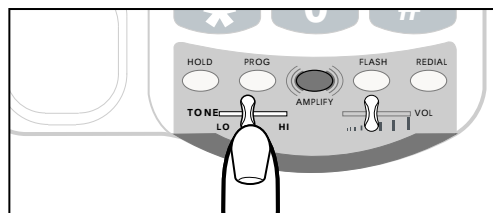
Figure 11 - Tone Slide Control

Using Your XL-30 (Continued)

Ringer Volume and Pitch

1. You can adjust the **RINGER VOLUME** as high as 85+ decibels. The settings available are **OFF**, **LO** and **HI** as shown in Figure 12.

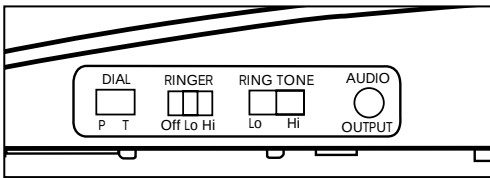
Figure 12 - Ringer Volume & Pitch

2. The **RING TONE** has two pitch settings: **LO** and **HI**.

Hold

To put the line on hold, press the **HOLD** button and return the handset to its cradle. You will notice the **HOLD** indicator light comes on.

NOTE: If you pick up an extension phone on the same line when the XL-30 is on **HOLD**, the XL-30 will release the **HOLD** and you can continue your conversation on the extension phone.

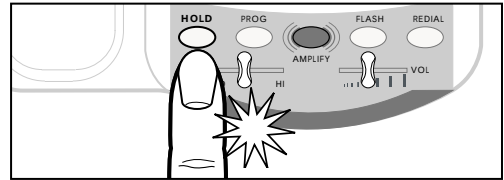
Figure 13 - Hold Button and Indicator

Last Number Redial

To redial the last number called, lift the handset and press **REDIAL**. The phone will redial up to 31 digits.

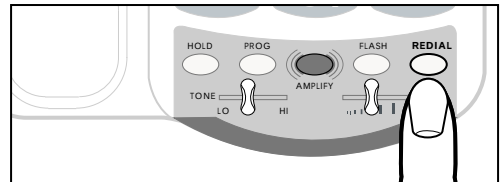
Figure 14 - Redial Button

Flash

Press the **FLASH** button shown in Figure 15 to access special services available from your local phone company such as Call Waiting and 3-Way Calling.

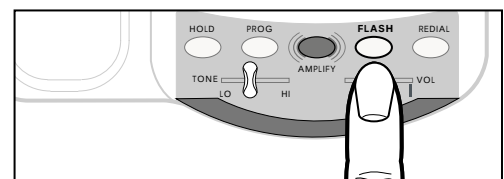
Figure 15 - Flash Button

Programming the Memory Buttons

You can automatically dial twelve programmed telephone numbers using the memory buttons shown in Figure 16. Each button can be programmed (or reprogrammed using the same procedure) for a phone number up to 15 digits long.

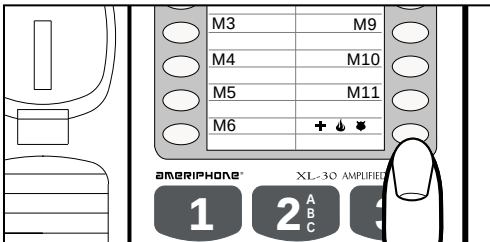
Figure 16 - Memory Buttons

To program a button, follow these steps:

1. Lift the handset and press **PROG** as shown in Figure 17.

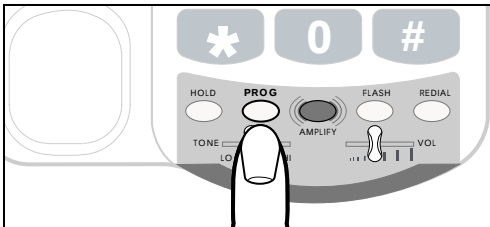
Figure 17 - PROG Button

2. Enter the phone number on the keypad (Figure 18) just as you would dial it normally.

Note: Press **REDIAL** if you want to insert a brief pause between two numbers.

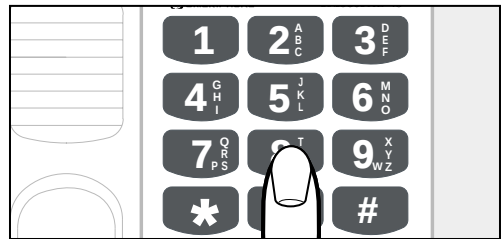
Figure 18 - Telephone Keypad

3. Press **PROG** (Figure 17) once again.
4. Press a memory button (Figure 16) to store this number on that button.
5. Immediately hang up the handset.
Note: Any number previously stored at that button number will be overwritten.
6. To change a stored number, repeat the programming process starting with step 1.
Note: Do not program 911 in any memory button.

Special Features (continued)

Direct Audio Output

The **AUDIO OUTPUT** socket allows you to connect the XL-30 to a hearing aid, neck loop, cochlear implant or other assistive listening devices.

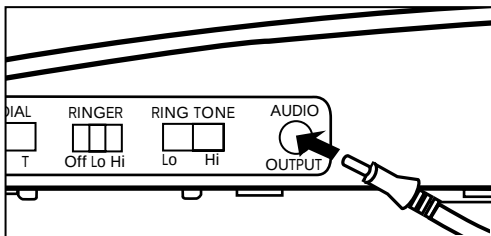
Figure 19 - Audio Output Socket

To use this port, plug in a compatible auxiliary cable and connect it to your assistive listening device. Speak into the handset when you use the Audio Output connection.

Outgoing Voice Amplification

The XL-30 can amplify your outgoing voice by up to 15 dB. To use this feature, set the outgoing voice amplifier switch to **AMPLIFY** (underneath the unit). **Note:** Even with the **AMPLIFY** switch **ON**, the voice of your callers can still be amplified by up to 15dB.

Note: The factory pre-setting is **NORMAL**.

Cochlear Implant Adapter Cord


A cochlear adapter cord is available from Ameriphone that can be used to plug the XL-30 directly into a speech processor.

WARNING! During thunderstorms, avoid using telephones, electric shock from lightning can occur. Before plugging any cord into your speech processor, you must consult your implant manufacturer.

Special Features (continued)

Connecting the Speech Processor:

- Plug the mono end of the plug into the speech processor.
- Plug the stereo end of the plug into the Audio Output socket on the phone.
- Speak into the handset when using the Cochlear Implant Adapter Cord.

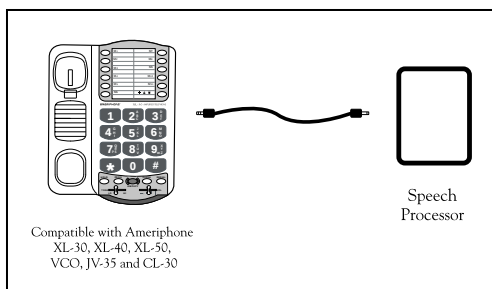Troubleshooting Chart

The chart below will help you solve most problems that may arise during operation of your telephone. Should the difficulty continue, contact Ameriphone or your authorized dealer for assistance.

SYMPTOM	CAUSE AND CORRECTIVE ACTION
No dial tone	Check all phone cord connections.
No number dialed when memory button pressed	No number stored in that memory button. See: "Programming the Memory buttons", P. 9.
Call cannot be dialed, or dials very slowly	See "Setting Up Your Phone" to reset the dial mode switch, P. 4, 5.
Call comes in but phone does not ring.	1. Check that line cord is not loose. 2. Check battery
Ring flasher not very bright	Check and replace the 9 Volt battery

Part 68 of FCC Rules Information

This equipment complies with Part 68 of the FCC rules and the requirements adopted by the ACTA. On the bottom of this equipment is a label that contains, among other information, a product identifier in the format US:AAAEQ##TXXXX. If requested, this number must be provided to the telephone company.

A plug and jack used to connect this equipment to the premises wiring and telephone network must comply with the applicable FCC Part 68 rules and requirements adopted by the ACTA. A compliant telephone cord and modular plug, RJ11 USOC, is provided with this product. It is designed to be connected to a compatible modular jack that is also compliant. See installation instructions for details.

The REN is used to determine the number of devices that may be connected to a telephone line. Excessive RENs on a telephone line may result in the devices not

ringing in response to an incoming call. In most but not all areas, the sum of RENs should not exceed five (5.0). To be certain of the number of devices that may be connected to a line, as determined by the total RENs, contact the local telephone company. For products approved after July 23, 2001, the REN for this product is part of the product identifier that has the format US:AAAEQ##TXXXX. The digits represented by ## are the REN without a decimal point (e.g., 03 is a REN of 0.3). For earlier products, the REN is separately shown on the label.

If this telephone equipment causes harm to the telephone network, the telephone company will notify you in advance that temporary discontinuance of service may be required. But if advance notice isn't practical, the telephone company will notify the customer as soon as possible. Also, you will be advised of your right to file a complaint with the FCC if you believe it is necessary.

The telephone company may make changes in its facilities, equipment, operations or procedures that could affect the operation of the equipment. If this happens the telephone company will provide advance notice in order for you to make necessary modifications to maintain uninterrupted service.

If trouble is experienced with this telephone equipment, for repair or warranty information, please contact Walker / Ameriphone, 1-800-874-3005. If the equipment is causing harm to the telephone network, the telephone company may request that you disconnect the equipment until the problem is resolved.

This telephone equipment is not intended to be repaired and it contains no repairable parts. Opening the equipment or any attempt to perform repairs will void the warranty. For service or repairs, call 1-800-874-3005.

Connection to party line service is subject to state tariffs. Contact the

state public utility commission, public service commission or corporation commission for information.

If your home has specially wired alarm equipment connected to the telephone line, ensure the installation of this telephone equipment does not disable your alarm equipment. If you have questions about what will disable alarm equipment, consult your telephone company or a qualified installer.

This telephone equipment is hearing aid compatible.

Customer-Owned Coin / Credit Card Phones:

To comply with state tariffs, the telephone company must be given notification prior to connection. In some states, the state public utility commission, public service commission or corporation commission must give prior approval of connection.

Part 15 of FCC Rules Information

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Your telephone equipment has been tested and found to comply with the limits of a Class B digital device, pursuant to Part 15 of FCC rules. These limits are designed to provide reasonable protection against harmful interference in residential installation.

This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation; if this

equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, you are encouraged to try to correct the interference by one of the following measures:

1. Where it can be done safely, reorient the receiving television or radio antenna.
2. To the extent possible, relocate the television, radio or other receiver with respect to the telephone equipment. (This increases the separation between the telephone equipment and the receiver.)
3. Connect the telephone equipment into an outlet on a circuit different from that to which the television, radio, or other receiver is connected.
4. Consult the dealer or an experienced radio/TV technician for help.

CAUTION: Changes or modifications not expressly approved by the manufacturer responsible for compliance could void the user's authority to operate the equipment.

Industry Canada Technical Specifications

This product meets the applicable Industry Canada technical specifications.

Before installing this equipment, users should ensure that it is permissible to be connected to the facilities of the local telecommunications company. The equipment must also be installed using an acceptable method of connection. In some cases, the company's inside wiring associated with a single line individual service may be extended by means of a certified connector assembly (telephone extension cord). The customer should be aware that compliance with the above conditions may not prevent degradation of service in some situations.

Repairs to certified equipment should be made by an authorized Canadian maintenance facility designated by the supplier. Any repairs or alterations made by the user to this equipment, or equipment malfunctions, may give the telecommunications company cause to request the user to disconnect the equipment.

Users should ensure for their own protection that the electrical ground connections of the power utility, telephone lines and internal metallic water pipe system, if present, are connected together. This precaution may be particularly important in rural areas.

Caution: Users should not attempt to make such connections themselves, but should contact the appropriate electric inspection authority, or electrician, as appropriate.

The Ringer Equivalence Number is an indication of the maximum number of terminals allowed to be connected to a telephone interface.

The termination on an interface may consist of any combination of devices subject only to the requirement that the sum of the Ringer Equivalence Numbers of all the devices does not exceed five.

[The term “IC:” before the certification / registration number only signifies that the Industry Canada technical specifications were met.]

SUPPLIER'S DECLARATION OF CONFORMITY

Issue Date: 09/23/02

We:

Walker, a Division of Plantronics
12082 Western Ave.
Garden Grove, CA 92841-2913

Declare that the product:

Description: Amplified Telephone

Models: XL-30

Bearing ACTA labeling identification number US:4B5TF09BXL-30
to which this declaration relates is in conformity with
the following standard/s:

Federal Communications Commission's 47 CFR Part 68

The Bases on which Conformity is being declared:

the Administrative Council on Terminal Attachments (ACTA)
adopted criteria:

TIA/EIA/IS-968 July 2001

Additional product information can be found on the website
www.ameriphone.com

A copy of this statement can be found at www.sdoc.ws

Acknowledgment of this declaration can be found on the ACTA
Website www.part68.org

Related Labeling Information:

AC REN: 0.9B

Jack Type: USOC type RJ11

This Phone is Hearing Aid Compatible

This mark first applied:

Name: LYNDALL FRY

Title: VP, Corporate Quality

Place of Issue:
Fremont Labs
42307 Osgood Rd Suite E
Fremont, CA USA 94539

Issued: 09/23/02

Warranty

This warranty applies only to Ameriphone products that are purchased and used in the United States or Canada.

Ameriphone warrants the XL-30 telephone against any defect in materials or workmanship for the period of one year from the date of purchase.

If your Ameriphone product is defective and returned within 30 days of the date of purchase, your Ameriphone dealer will replace it at no charge.

If returned after 30 days but within one year from the date of purchase, we will repair or replace it at no charge. In the repair of your XL-30, we may use new or reconditioned replacement parts. If we elect to replace your XL-30, we may replace it with a new or reconditioned product of the same or similar design. Repair or replacement will be warranted for either 90 days or the remaining time on the original warranty period, whichever is longer.

Implied warranties, including those of fitness for a particular purpose and merchantability (an unwritten warranty that the product is fit for ordinary use), are limited to one year from date of purchase. We will not pay for loss of time, inconvenience, loss of use of your XL-30 phone, or property damage caused by your XL-30 phone or its failure to work, or any other incidental or consequential damages. Some states do not allow limitations on how long an implied warranty lasts or the exclusion of incidental or consequential damages, so the above exclusions or limitations may not apply to you.

To get warranty service for your XL-30 phone, you must provide proof of the purchase date. Within 30 days of the date of purchase, return your XL-30 to the place where you purchased it for immediate replacement. After 30 days, call Ameriphone at (800) 874-3005 voice or 800-772-2889 TTY for the authorized service center nearest you. You

Warranty

must prepay all shipping costs. We suggest you save the original package materials in the event you need to ship the XL-30. When shipping for warranty repair, include your name, address, phone number, proof of date of purchase, and a description of the problem. After repairing the product, we, (or the service center) will ship it back to you at no cost within the United States and Canada. CANADIAN RESIDENTS: call ameriphone at 800-874-3005, 800-772-2889 TTY for instructions.

This warranty does not cover defects resulting from accidents, damage while in transit to our service location, alterations, unauthorized repair, failure to follow instructions, misuse, use outside the United States or Canada, fire, flood, and acts of God. Nor do we warrant the product to be compatible with any particular telephone equipment, party line, key telephone systems or more sophisticated switching systems. If your XL-30 is not

covered by this warranty, call us at (800) 874-3005 voice or (800) 772-2889 TTY for advice as to whether we will repair your XL-30 and other repair information. The repair shall be warranted for 90 days.

Specifications

Maximum gain

High frequency (3kHz): 30dB

Wide band (300 to 3000Hz):
20dB

Tone control range

Full range: 300 to 3000Hz using
one slide control.

Dimensions

Size: 9 ¹/₂" x 7" x 3 ¹/₄"

Weight: 2.52 lbs.

Power Requirements

One 9 Volt alkaline battery for
ring flasher operation only. (not
included)

Connecting for Desk or Wall Mount	
Desktop Use	4
Wall Mount	5
Direct Audio Output	10
Features	3
Flash	8
Help from Ameriphone	1
Hold	8
Incoming Voice Tone	7
Incoming Voice Volume	6
Installing the Backup Batteries ..	4
Last Number Redial	8
Programming the	
Memory Buttons	9
Ringer Volume and Pitch	8
Sales Receipt	1
Troubleshooting	12
Warranty Service	1

Ameriphone products are manufactured by Walker, a division of Plantronics, Inc.
© 2002 Plantronics, Inc. All rights reserved. Ameriphone, Amplifying Your Life, and Walker are trademarks or registered trademarks of Plantronics, Inc.

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